Bioclin

AMILASE CINÉTICA

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FINALIDADE

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- Winn-Deen ES, David H, Sigler G and Chavez R. Development of a direct assay for a-amylase. Clin Chem 1988; 34: 2005-2008.
- Balsells D, Gella FJ, Gubern G, Canalias F. Reference values for a-amylase in human serum and urine using 2-Chloro-4-Nitrophenyl-α-D-Maltotrioside as substrate. Clin Chim Acta 1998; 274: 213-217.
- Tietz Textbook of Clinical Chemistry, 2nd edition. Burtis CA, Ashwood ER. WB Saunders Co., 1994.
- WHO. Use of anticoagulants in diagnostic laboratory investigations and stability of blood, plasma and serum samples. WHO/DIL/LAB/99.1 rev. 2, 2002:31.
- QUIBASA: Dados do Departamento de Pesquisa e Desenvolvimento.

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INSTRUCTIONS FOR USE


FUNCTION

Method for determination of the α -Amylase in biological samples of serum, heparinized plasma, urine and cavity fluids (ascitic, duodenal or pleural). Kinetic test, only for *in vitro* diagnostic use.

PRINCIPLE OF ACTION

Methodology: Kinetic CNPG

The α -Amylase catalyzes the hydrolysis of the α -(2-Chloro-4-Nitrophenyl)- β -1,4-Galactopiranosilmaltoside (Gal-G2- α -CNP) releasing 2-Chloro-4-Nitrophenol (CNP) and 1,4 Galactopiranosilmaltoside (Gal-G2). The catalyst concentration is determined from the rate of 2-Chloro-4-Nitrophenol formation, measured at 405 nm.

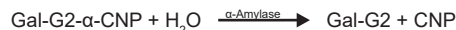REAGENTS

Number 1 - Buffered Substrate - Store between 2 and 8°C. Contains: MES buffer < 200 mmol/L, Gal-G2- α -CNP < 5 mmol/L, activators, stabilizers and preservative.

PRESENTATION

Presentation	Reagent N° 1
1	3 x 20 mL
2	5 x 20 mL
3	5 x 40 mL

EQUIPMENTS AND OPERATIONAL INPUTS

Thermostated spectrophotometer, pipettes, watch or stopwatches, test tubes, Biocontrol N, Biocontrol P Bioclin and Biocal Bioclin. They can be found at markets specialized on Laboratories of Clinical Analysis.

TRANSPORTATION AND STORAGE CONDITIONS

The storage temperature should be between 2 to 8°C. The transport at temperatures up to 30°C should not exceed 5 days. Protect from light and avoid moisture. **Do not freeze.**

SPECIAL CARE

1- For *in vitro* diagnostic use only.

2- Strictly follow the methodology proposed to obtain exact results.

3- Water used in material cleaning must be recent and free of contaminants.

4- Saturated deionizer columns release alkaline water, many ions, oxidizing agents and reducers that may alter the results significantly.

5- Visible hemolysis may cause slight variations in results.

6- Reagent N°1 is easily contaminated by saliva. Therefore, we recommend the use of masks, pears and automatic pipettes for handling this product.

7- Reagent N° 1 contain Sodium Azide which can be irritating to skin and mucous membrans. Handle with care.

8- We recommend applying the local, state and federal rules for environmental protection, so that disposal of reagents and biological material can be made in accordance with current legislation.

9- To obtain information related to biosafety or in case of accidents with the product, consult the SDS (Safety Data Sheet) available on the website www.bioclin.com.br or upon request by the SAC (Customer Advisory Service) of Quibasa.

10- Do not use the product in case of damaged packaging.

11- It is essential that the instruments and equipments used are properly calibrated and subjected to periodic maintenance.

SAMPLES

Serum obtained free of hemolysis, plasma with Heparin, Urine and cavity fluids (ascitic, duodenal and pleural). Other anticoagulants (EDTA, Citrate, Oxalate) inhibit the activity of Amylase.

Serum or Plasma: stable for 7 days at 15 to 30°C, 2 months at 2 to 8°C and 1 year at -20°C.

Urine: stable for 2 days at 15 to 30°C, 10 days at 2 to 8°C and 3 weeks at -20°C. Gather in the range of 2 to 24 hours.

Cavity fluids (ascitic, duodenal and pleural): Stable 7 days at -20°C.

PROCESS DESCRIPTION

The stability of the Amylase Kinetic product after opening (In Use / On Board – Product installed in equipment) is 28 days. This stability may vary according to test conditions. Therefore, it is suggested to monitor the performance of the product using control serum.

REAGENTS PREPARATION

The Working Reagent is ready for use.

REACTION CONDITIONS

It is indispensable the use of thermostated cuvette at 37°C, light path of 1cm and reading at 405 nm.

TECHNIQUE

Bioclin recommends using the kit Biocal as calibrator and as control serum Biocontrol N and P Bioclin kits.

Serum or Plasma

Add 20 μ L of sample to 1.0 mL of Buffered Substrate, mix and transfer to a thermostated cuvette at 37°C. Make the initial reading, simultaneously starting the timer. Repeat the readings at 1, 2 and 3 minutes. Calculate the mean difference in absorbance per minute ($\Delta A/\text{min.}$) and used to calculate the result.

Urine

Mix, measure and note the urinary volume in mL and the time of sample collection. Proceed dosage as described for serum.

CALCULATIONS

Serum or Plasma

Amylase (U/L) = $\Delta A/\text{min.} \times 3954$

Results are expressed in U/L.

Urine

Amylase (U/L) = $\Delta A/\text{min.} \times 3954$

Urinary Amylase (U/h) = $\frac{\text{Amylase (U/L)} \times V}{(H \times 1000)}$

V = Urinary volume in mL

H = Hours from urine collection

Example

Urinary Amylase (value found in the dosage) = 90 U/L

V = 200 mL

H = 2

Urinary Amylase (U/h) = $(90 \times 200) / (2 \times 1000) = 9 \text{ U/h}$

Amylase clearance / creatinine clearance ratio

Evaluation of the amylase clearance / creatinine clearance ratio is important for determining the diagnosis of acute pancreatitis. To define the clearance ratio (Amylase / Creatinine) of the patient, amylase activity and creatinine concentration in serum and a urine sample should be determined.

Equation:

$$\text{Ratio (\%)} = \frac{\text{Urine Amylase (U/L)} \times \text{Serum Creatinine (mg/dL)}}{\text{Serum Amylase (U/L)} \times \text{Urine Creatine (mg/dL)}} \times 100$$

For an average variation in absorbance at 405 nm greater than 0.25, repeat the determination, diluting the sample with 0.85% NaCl. Multiply the result by the dilution factor. To convert the data in the International System (SI), multiply the result: U/L x 0.01667 = $\mu\text{kat/L}$

PROCEDURE LIMITATIONS

The specifications below are for semi-automatic equipment: The kinetic method is based on the molar absorptivity, therefore, the readings should be performed in a spectrophotometer that satisfies the following conditions: Wavelength 405 nm

Semi trajectory of the pass band 10 nm.

Stray light less than 0.5%.

1cm thermostated cuvette

INTERFERING

Saliva contains amylase. Therefore, you should not pipet by mouth and avoid contact of the reagent with the skin. The lipemia (triglycerides up to 1000 mg/dL) and bilirubin (up to 20 mg/dL) do not interfere with the methodology. Hemoglobin up to 2.5 g/L no affect the results. Some drugs may interfere with laboratory results, elevate the levels of serum amylase (morphine, meperidine codeine, thiazide diuretics) or decreasing, as in cases of barbiturate poisoning.

INTERNAL QUALITY CONTROL

The Clinical Laboratory must have an internal quality control, where all procedures, rules, limits and tolerance to variations be clearly established. It is important to mention that all measurement systems present a analytical variety, and it must be monitor by the laboratory. Therefore, it is recommendable the use of controls, allowing the precision and accuracy of the dosages.

TRACEABILITY

The kit calibration can be made using the theoretical calibration factor based on the molar absorptivity of the 2-Chloro-4-Nitrofenol, or through the BIOCAL calibrator. Bioclin recommends the usage of the BIOCAL calibrator which is traceable to the reference material IRMM/IFCC-456.

REFERENCE VALUES

The reference values in U/L for this method were obtained by determining the α-Amylase in healthy populations of males and females in various age groups.
Serum: < 125 U/L
Urine: < 450 U/L or < 30 U/h

Amylase clearance / creatinine clearance ratio: 1.0 to 4.0 %

These values should be used as guidelines, each laboratory should establish its range of reference values, according to the population served.
The results provided by this kit must be interpreted by the medical professional responsible, not being the only criterion for the determination of diagnosis and/or treatment of the patient.

PRODUCT PERFORMANCE

ACCURACY

Comparison of Methods

The Amylase Kinetic kit was compared with another method commercially available to measure Amylase. 42 analysis were conducted and the results were evaluated. The linear equation was obtained $Y = 1.0771X + 0.3224$, with a correlation coefficient equal to 0.9998. With these results we can conclude that the kit has good methodological specificity.

PRECISION

Repeatability

The repeatability was calculated from 40 successive determinations, using 3 samples with different concentrations, obtaining the following results:

Repeatability	Sample 1	Sample 2	Sample 3
Average Concentration (U/L)	110.33	277.60	194.63
Standard Deviation (U/L)	1.05	2.53	1.44
Coefficient of Variation (%)	0.95	0.91	0.74

Reproducibility

The reproducibility was calculated from 40 successive determinations for 3 consecutive days, using 3 samples with different concentrations, obtaining the following results:

Reproducibility	Sample 1	Sample 2	Sample 3
Average Concentration (U/L)	111.73	276.70	194.52
Standard Deviation (U/L)	1.23	3.46	2.56
Coefficient of Variation (%)	1.10	1.25	1.32

SENSITIVITY

The sensitivity was calculated from 40 determinations of a sample free of Amylase. The average found was of 2.23 U/L with standard deviation of 0.42 U/L. The sensitivity, which indicates the method detection limit, corresponds the average plus 3 times the standard deviation and is equal to 3.49 U/L.

LINEARITY

The reaction is linear up to 2000 U/L.

DIAGNOSTIC SIGNIFICANCE

The Amylase, predominantly of pancreatic and salivary, is normally present in blood and urine in small quantities. Increases rapidly in plasma after the onset of acute pancreatitis, where there is a parallel increase in urinary amylase.
Elevations are also found in mesenteric infarction, perforated gastric ulcer, carcinoma of the head of the pancreas, mumps, kidney failure, diabetic acidosis. Decreased plasma levels are observed in hepatitis, liver cirrhosis, toxemia of pregnancy, eclampsia, pancreatic carcinoma.

BIBLIOGRAPHIC REFERENCES

- Winn-Deen ES, David H, Sigler G and Chavez R. Development of a direct assay for a-amylase. Clin Chem 1988; 34: 2005-2008.
- Balsells D, Gella FJ, Gubern G, Canalias F. Reference values for a-amylase in human serum and urine using 2-Chloro-4-Nitrophenyl-a-D-Maltotrioside as substrate. Clin Chim Acta 1998; 274: 213-217.
- Tietz Textbook of Clinical Chemistry, 2nd edition. Burtis CA, Ashwood ER. WB Saunders Co., 1994.
- WHO. Use of anticoagulants in diagnostic laboratory investigations and stability of blood, plasma and serum samples. WHO/DIL/LAB/99.1 rev. 2, 2002:31.
- QUIBASA: Dados do Departamento de Pesquisa e Desenvolvimento.

QUALITY ASSURANCE

Before being released for consumption, all Bioclin reagents are tested by the Department of Quality Control. The quality of reagents is assured until expiration date stated on the presentation packaging, when stored and transported under appropriate conditions.

QUIBASA QUÍMICA BÁSICA Ltda
Rua Teles de Menezes, 92 - Santa Branca
CEP 31565-130 - Belo Horizonte/MG - Brasil
Phone.: +55 31 3439.5454 | E-mail: bioclin@bioclin.com.br
CNPJ: 19.400.787/0001-07 - Made in Brazil

CUSTOMER SERVICE
Customer Advisory Service
Phone.: 0800 0315454 | E-mail: sac@bioclin.com.br

ANVISA registration for Amylase Kinetic kit: 10269360079

Review: August/2024

UNIVERSAL SYMBOLOGY

CATALOG NUMBER

LOT NUMBER

MANUFACTURING DATE

VALIDITY DATE
(last day of the month)

TEMPERATURE LIMIT
(store)

CONTENT IS SUFFICIENT
FOR <N> TEST

SEE INSTRUCTIONS
FOR USE

IN VITRO DIAGNOSTIC
PRODUCT

KEEP AWAY
FROM SUNLIGHT

DO NOT REUSE

CAUTION

MADE BY

CONTROL

CONTROL +
POSITIVE CONTROL

CONTROL -
NEGATIVE CONTROL

BIOLOGICAL RISK

FLAMMABLE

CORROSIVE

TOXIC

DO NOT USE IF
PACKAGE IS
DAMAGED

STERILE
PRODUCT
STERILIZED

DANGER