CK NAC UV

REF K010

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REF K010

INSTRUCTIONS FOR USE

TO OBTAIN THE INSTRUCTIONS FOR USE IN PRINTED FORMAT, AT NO ADDITIONAL COST, CONTACT CUSTOMER ADVISORY SERVICE:

SAC: +55 (31) 3439 5454 / 0800 031 5454 / sac@bioclin.com.br


FUNCTION

Method for determination of the creatine kinase (CK) in biological samples of serum or plasma (EDTA or Heparin). Kinetic test, only for *in vitro* diagnostic use.

PRINCIPLE OF ACTION

Methodology: Kinetic UV

Determination of creatine kinase according to the reactions:

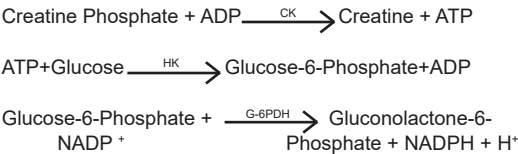CK - Creatine Kinase
G-6-PDH – Glucose-6-Phosphate-Dehydrogenase
HK - Hexokinase

The reduction speed of NADP⁺ to NADPH is proportional to the activity of the CK in the sample.

REAGENTS

Number 1 - Buffer - Store between 2 and 8°C. Contains: Imidazole Buffer < 300 mmol/L, Glucose < 100 mmol/L, NADP < 10 mmol/L, Hexokinase < 50 KU/L, N-Acetylcysteine < 100 mmol/L, stabilizers, surfactant and preservative.

Number 2 - Enzyme Substrate - Store between 2 and 8°C. Contains: Buffer < 300 mmol/L, Glucose-6-Phosphate-Dehydrogenase < 50 KU/L, Creatine Phosphate < 300 mmol/L, ADP < 30 mmol/L, stabilizers, surfactant and preservative.

PRESENTATION

Reagent	Presentation				
	1	2	3	4	5
Nº 1	1 x 40 mL	2 x 40 mL	4 x 40 mL	2 x 40 mL	4 x 40 mL
Nº 2	1 x 10 mL	2 x 10 mL	4 x 10 mL	1 x 20 mL	2 x 20 mL

EQUIPMENTS AND OPERATIONAL INPUTS

Thermostated spectrophotometer, pipettes, test tubes, watch or stopwatch, Biocal Bioclin, Biocontrol N and Biocontrol P Bioclin. These itens are found at markets specialized on Laboratories of Clinical Analysis.

TRANSPORTATION AND STORAGE CONDITIONS

The storage temperature should be between 2 to 8°C. The transport at temperature up to 30°C should not exceed 5 days. Protect from light and avoid moisture. **Do not freeze.**

SPECIAL CARE

- 1- For *in vitro* diagnostic use only.**
- Strictly follow the methodology proposed to obtain exact results.
- Water used in material cleaning must to be recent and free of contaminants.
- Saturated deionizer columns release alkaline water, many ions, oxidizing agents and reducers that may alter the results significantly.
- 5- It is important for the proper performance of the test, a strict control of time, temperature and pH.**
- We recommend applying the local, state and federal rules for environmental protection, so that disposal of reagents and biological material can be made in accordance with current legislation.
- To obtain information related to biosafety or in case of accidents with the product, consult the SDS (Safety Data Sheet) available on the website www.bioclin.com.br or upon request by the SAC (Customer Advisory Service) of Quibasa.
- Do not use the product in case of damaged packaging.
- It is essential that the instruments and equipments used are properly calibrated and subjected to periodic maintenance.

SAMPLES

Serum or plasma collected with Heparin or EDTA. The sample is stable for 4 hours between 15 and 30°C ³, 7 days between 2 and 8°C and 30 days at -20°C ³. Do not use hemolyzed sample or lipemic.

PROCESS DESCRIPTION

The calibration stability of the CK NAC UV kit installed on refrigerated equipment is at least 11 days. This stability may vary depending on the conditions of the test, equipment and environment. Therefore, it is suggested to follow the product performance using control serum.

PREPARATION OF WORKING REAGENT

Mix 4 parts of reagent Nº 1 with 1 part of reagent Nº 2. Stable for 2 days at temperatures between 15 to 30°C and 15 days if kept between 2 to 8°C. Do not use hemolyzed sample.

REACTION CONDITIONS

It is indispensable the use of a thermostated cuvette at 37°C, light path of 1cm and reading at 340 nm.

TECHNIQUE

Bioclin recommends using the kit Biocal Bioclin as calibrator and as control serum, Biocontrol N and P Bioclin Kits. To do this, mix the reagents according to the information in the table below.

Working Reagent	1.0 mL
Sample	20 µL

Mix and transfer to a thermostated cuvette at 37°C and wait 2 minutes. Make the initial reading, triggering the timer. Repeat the readings at 1, 2 and 3 minutes.

CALCULATIONS

Calculate the mean difference in absorbance per minute (Δ A/min.) and use this value to calculate the result.

CK (U/L) 340 nm = Δ A/min x 8095

PROCEDURE LIMITATIONS

The readings should be taken in a spectrophotometer that satisfies the following conditions:
Wavelength 340 nm
Semi trajectory of the pass band 10 nm
Spurious light lower than 0.5%
10 mm Thermostated cuvette.

INTERFERENT

No interference was observed by Glucose until 7 g/L, Hemoglobin until 5 g/L and Triglycerides until 7 mmol/L.

INTERNAL QUALITY CONTROL

The Clinical Laboratory must have an internal quality control, where all procedures, rules, limits and tolerance to variations be clearly established. It is important to mention that all measurement systems present a analytical variety, and it must be monitor by the laboratory. Therefore, it is recommendable the use of controls, allowing the precision and accuracy of the dosages.

TRACEABILITY

The kit calibration can be made using the theoretical calibration factor based on the molar absorptivity of the NADP, or through the BIOCAL calibrator which is traceable to the reference material ERM-AD455 / IFCC and the IFCC reference method.

REFERENCE VALUES

The reference values for this method were obtained by determining CK in healthy populations of male and female.

Adults (37°C)	
Male	24 - 195 U/L
Female	24 - 170 U/L

Children and Teenagers (37°C)		
Age	Male	Female
1 - 3 month	29 - 303 U/L	40 – 474 U/L
4 - 12 month	25 – 172 U/L	27 – 242 U/L
13 - 24 month	28 – 162 U/L	25 – 177 U/L
2 - 10 years	31 – 152 U/L	25 – 177 U/L
11 - 14 years	31 – 152 U/L	31 – 172 U/L
15 - 18 years	34 – 147 U/L	28 – 142 U/L

These values should be used as guidelines, each laboratory should establish its range of reference values, according to the population served.

The results provided by this kit must be interpreted by the medical professional responsible, not being the only criterion for the determination of diagnosis and/or treatment of the patient.

PRODUCT PERFORMANCE
ACCURACY

Comparison of Methods

The CK NAC UV kit was compared to other method for CK NAC dosage available in the market. 42 analysis were performed and the results were evaluated. The linear equation obtained was $Y = 1.0045x + 2.5661$, with linear correlation coefficient equal to 0.9994. With these results we can conclude that the kit shows good methodological specificity.

PRECISION

Repeatability

The repeatability was calculated from 40 successive determinations, using 3 samples with different concentrations, obtaining the following results:

Repeatability	Sample 1	Sample 2	Sample 3
Average Concentration(U/L)	72.55	412.28	235.85
Standard Deviation (U/L)	0.88	2.89	1.64
Variation Coefficient (%)	1.21	0.70	0.70

Reproducibility

The reproducibility was calculated from 40 successive determinations for 3 consecutive days, using 3 samples with different concentrations, obtaining the following results:

Reproducibility	Sample 1	Sample 2	Sample 3
Average Concentration(U/L)	72.65	411.85	234.68
Standard Deviation (U/L)	0.98	1.61	5.24
Variation Coefficient (%)	1.35	0.39	2.23

SENSITIVITY

The sensitivity was calculated from 40 determinations of a sample free of the analyte (CK NAC). The mean concentration of CK NAC UV was 0.225 U/L with standard deviation of 0.480 U/L. The sensitivity, which indicates the limit of detection, which corresponds to the average plus three times the standard deviation is equal to 1.664 U/L.

LINEARITY

The method is linear to 2000 U/L. If the variation of molar extinction per minute is larger than 0,250, repeat the test using a dilution of 1:10 with physiologic serum. The result obtained is multiplied by 10.

DIAGNOSTIC SIGNIFICANCE

The enzyme Creatine Phosphokinase (CPK) and Creatine Kinase (CK) is found in high concentrations in cardiac and skeletal muscles and brain.

The values of CK are elevated in the presence of muscle damage as in myocardial infarction. Following the attack, the serum concentration of this enzyme increases rapidly (3 to 6 hours), reaching a peak after 12 to 24 hours and remains elevated for a short period of 2 to 3 days. Other potential sources of elevated total CK are: hypothyroidism, muscle disease (myopathy, polymyositis), stroke, seizures, cardiac catheterization, multiple trauma, strenuous and prolonged immobilization.

BIBLIOGRAPHIC REFERENCES

- 1 -SCANDINAVIAN Society For Clinical Chemistry Scand. J. Clin. Lab. INVEST., 1974, 33, 291.
- 2 - Ann. Biol. Clin., 1982, 40,99.
- 3 - WHO. Use of anticoagulants in diagnostic laboratory investigations and stability of blood, plasma and serum samples. WHO/DIL/LAB/99.1 rev. 2, 2002:31
- 4 - QUIBASA: Dados do Departamento de Pesquisa e Desenvolvimento.

QUALITY ASSURANCE

Before being released for consumption, all **Bioclin** reagents are tested by the Department of Quality Control. The quality of reagents is assured until expiration date stated on the presentation packaging, when stored and transported under the appropriate conditions.

QUIBASA QUÍMICA BÁSICA Ltda

Rua Teles de Menezes, 92 - Santa Branca
CEP 31565-130 - Belo Horizonte - MG - Brasil
Phone.: +55 31 3439.5454 | E-mail: bioclin@bioclin.com.br
CNPJ: 19.400.787/0001-07 - Made in Brazil

CUSTOMER SERVICE

Customer Advisory Service
Phone.: 0800 0315454 | E-mail: sac@bioclin.com.br

ANVISA registration for CK NAC UV kit: 10269360092

Review: October/2024

UNIVERSAL SYMBOLOGY

CATALOG NUMBER

LOT NUMBER

MANUFACTURING DATE

VALIDITY DATE
(last day of the month)

TEMPERATURE LIMIT
(store)

CONTENT IS SUFFICIENT
FOR <N> TEST

SEE INSTRUCTIONS
FOR USE

IN VITRO DIAGNOSTIC
PRODUCT

KEEP AWAY
FROM SUNLIGHT

DO NOT REUSE

CAUTION

MADE BY

CONTROL

POSITIVE CONTROL

NEGATIVE CONTROL

BIOLOGICAL RISK

FLAMMABLE

CORROSIVE

TOXIC

DO NOT USE IF
PACKAGE IS
DAMAGED

PRODUCT
STERILIZED

DANGER