Bioclin

ERLIQUIOSE IgG VET FAST

REF VET001

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FINALIDADE

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PURPOSE

Fast immunochromatographic test for the qualitative determination of anti-*Ehrlichia canis* IgG antibodies in serum, plasma or whole blood samples.

PRINCIPLE OF ACTION

Methodology: Immunochromatography.
Ehrlichiosis IgG VET FAST is an *in vitro* diagnostic test for the rapid and qualitative detection of anti-*Ehrlichia canis* IgG antibodies in blood, serum or plasma samples from dogs.
The cassette consists of a nitrocellulose membrane impregnated with the *Ehrlichia canis* recombinant antigen test line. The conjugate consists of a mouse antibody anti-dog IgG conjugated to colloidal gold particles. When the sample is added to the cassette, it migrates by capillary action until it reaches the conjugate. In this region, the IgG antibodies present in the sample interact with the conjugate, forming the complex: dog IgG – conjugate.
The sample continues to migrate through the immunochromatographic strip until it reaches the test line, where the dog IgG antibody also interacts with the immobilized antigen of *Ehrlichia canis* in the test line, forming the complex: *Ehrlichia canis* antigen – dog IgG – conjugate. This shows the formation of a pink-red line on the test line; indicating positive result.
Sample excess continues to migrate through the test strip until it reaches the control line, which is composed of impregnated anti-conjugate antibody. Thus, the formation of a second line (control line) is observed; indicating that the test yielded a valid result.

REAGENTS

Reagent Number 1 (R1): Cassette - Store between 15 and 30°C. **Do not freeze.** It contains:

- A)** 01 plastic cassette with 01 test strip containing
- Control region (C): Goat antibody anti-mouse antibody immobilized.
 - Test region (T): Immobilized *Ehrlichia canis* recombinant antigen.
 - Conjugate: Mouse antibody anti-dog IgG, colloidal gold conjugate.

B) 01 sachet of silica

Reagent Number 2 (R2): Diluent - Store between 15 and 30°C. **Do not freeze.** Contains: Buffer solution and preservative.

Collection Accessories:

- Disposable plastic pipette

PRESENTATIONS

Reagent	Number of Cassettes for Package	Diluent
Cassette + Diluent	10	1 x 2,5 mL
	25	1 x 5,0 mL
	10*	10 x 0,5 mL
	25*	25 x 0,5 mL
	5	1 x 1,25 mL

*Individual presentations.

EQUIPMENTS AND OPERATIONAL INPUTS

Clock or stopwatch. They are in the specialized market of articles for Clinical Analysis Laboratories.

TRANSPORTATION AND STORAGE CONDITIONS

The storage and transport temperature of the product must be between 15 and 30°C. Keep away from direct sunlight. **Do not freeze.** The cassette is sensitive to air humidity. Therefore, it should always be kept in its inviolate primary packaging (aluminized sachet) until the time of use.

SPECIAL CARE


- 1- Product for in vitro diagnosis, for veterinary use only.**
- 2- Strictly follow the methodology proposed for obtaining results.
- 3- Cassettes (test devices) do not contain contaminants; but since the samples are biological, the professional must strictly follow the safety rules and routines. The use of disposable gloves and other personal protective equipment is essential.
- 4- Use clean, dry containers to package biological samples.
- 5- Open the aluminum envelop to remove the cassette only at the time of use.
- 6- The kit cannot be frozen or exposed to high temperatures, since this will cause deterioration of its components.
- 7- The reagents, as well as the samples, should reach room temperature before starting the test.
- 8- Do not use reagents from different lots or kits.
- 9- Do not use kit beyond its expiration date. The expiry date is stated on the product label.
- 10- The cassette (test device) and the capillary tube are disposable and cannot be reused.
- 11- The kit materials that had contact with the sample of the patient must have a special destination. The capillary tube the cassettes (test devices) and diluent vials should be disposed of in organic waste.
- 12- Secondary packaging (box), Instruction for use, aluminum envelope and silica sachet can be disposed of in trash common.
- 13- We recommend to apply the local, state and federal environmental protection so that the disposal of reagents and biological material is done in accordance with the current legislation.
- 14- To obtain information related to biosafety or in case of accidents with the product, consult the MSDS (Material Safety Data Sheet) available on the website www.bioclin.com.br or through request by the Customer Service Department (SAC) of the Quibasa.

SAMPLES

Perform the test with serum, plasma or whole blood.
Serum: Allow blood to clot and then centrifuge to separate serum.
Plasma: Collect whole blood in a tube containing anticoagulants such as heparin, citrate, or EDTA. Centrifuge the blood and separate the plasma.
Whole blood: Collect the blood sample in a tube containing EDTA, heparin, or citrate using standard laboratory procedures.
Whole blood samples should be analyzed within 24 hours, can be stored at 2-8°C and never be frozen.
Serum and plasma samples should be analyzed immediately after collection, but if necessary, they can be stored at 2-8°C for up to 7 days or at -20°C for longer periods.

DESCRIPTION OF PROCESS

- 1- The sample should be at room temperature before starting the test.
 - 2- Remove the cassette (test device) from the protective packaging and place it on a clean, dry and level surface
 - 3- Collect the sample using the disposable capillary tube that accompanies the kit.
 - 4- Tighten the capillary tube above the marked trace and touch the sample. Relieve the pressure so the sample is aspirated. Aspirate until the trace marked to obtain the volume of 10 µL.
 - 5- Dispense the sample by pressing the pipette again, touching the test membrane at the sample well.
- Note: Alternatively, an automatic pipette can be used to collect and add the sample.
- 6- Hold the Diluent bottle (Reagent N° 2) vertically and apply 2 drops (70 µL) in sample well.
 - 7- Wait for the lines to form. Interpret the results between 5 and 10 minutes. Do not interpret after 10 minutes.

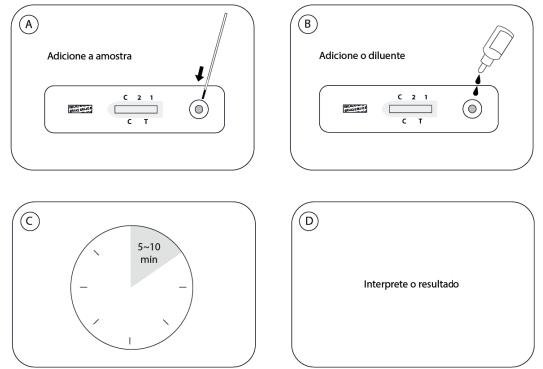INTERPRETATION OF RESULTS

Reactive Test: Formation of two lines: one in the control region (C) and one in the test region (T). Do not interpret after 10 minutes.	
Not Reactive Test: Line formation only in control region (C). Do not interpret after 10 minutes.	
Invalid Test: The absence of line formation in the control region (C), with or without line formation in the test region (T), indicates procedure error or deterioration of the cassette. In this case, repeat the test using a new cassette.	

TEST LIMITATIONS

This product is intended for screening for canine Ehrlichiosis. Interpretation should not be based on a single essay. Other confirmatory tests must be included before a sample is considered positive. Therefore, for the diagnosis of the disease, the professional must take into account not only the result of this test; but also, clinical signs, animal health conditions, and other diagnostic results that may be available.

INTERNAL QUALITY CONTROL

The Clinical Laboratory must have an internal program of quality control, where procedures, standards, limits and tolerance for variations are clearly established. It is important to note that all measuring systems characteristic analytical variability, which should be monitored by the laboratories themselves. Therefore, it is the use of controls, which allow the assessment of the precision and accuracy of dosages.

PRODUCT PERFORMANCE
QUALITY CONTROL

CLINICAL SENSITIVITY AND SPECIFICITY

This product has been tested comparing with other methods. A total of 95 samples characterized by canine Ehrlichiosis were analyzed; 37 samples were diagnosed as positive and 58 were diagnosed as negative. The results show that the clinical sensitivity of the product is 94,59% and the clinical specificity is 94,82%.

METHOD		REFERENCE		Total
Cassette Ehrlichiosis IgG VET FAST	Result	Reactive	No reactive	
	Reactive	35	3	
	Not Reactive	2	55	57
Total result		37	58	95

Clinical sensitivity: (35/37): 94,59%
Clinical Specificity: (55/58): 94,82%
Precision: (90/95): 94,73%

Precision
REPEATABILITY

Repeatability was calculated from 10 successive determinations, using 2 different samples, obtaining the following results:

Samples	N° of Repetitions	Expected Result	Found Result
Sample 1	10	Reactive	Reactive
Sample 2	10	Reactive	Reactive

REPRODUCIBILITY

Reproducibility was calculated from 10 successive determinations over 3 consecutive days, using 2 different samples, obtaining the following results:

Samples	N° of Days	N° de Repetitions	Expect- ed Result	Found Result
Sample 1	03	10	Reactive	Reactive
Sample 2	03	10	Reactive	Reactive

CLINICAL SIGNIFICANCE

Ehrlichiosis is a disease caused by bacteria of the *Ehrlichia* genus, which parasitize blood cells such as leukocytes, erythrocytes, endothelial cells and platelets. Transmission occurs through the bite of the brown tick (*Rhipicephalus sanguineus*). The tick is contaminated with the bacteria after the blood meal in dogs with high parasitemia, which commonly occurs in the acute phase of the disease. The incubation period for ehrlichiosis in dogs varies between 8 to 20 days, and the cynical signs observed include petechiae in the skin and mucous membranes, lymph node enlargement, hepatomegaly and splenomegaly, apathy, weight loss, inappetence, anorexia, fever, runny nose and eye. When treatment does not result in total elimination from the body, further relapses may occur. Rapid diagnosis based on the detection of antibodies against *Ehrlichia* allows the establishment of an adequate treatment protocol, avoiding the worsening of the disease.condition.

BIBLIOGRAPHIC REFERENCES

- 1- COUTO, C.G. Doenças Rickettsiais In: BIRCHARD, SHERDING, Manual Saunders: Clínica de pequenos animais. Ed. Roca: 139-42, 1998
2- ETTINGER, S. J. FELDMAN, E. C. Tratado de Medicina Interna Veterinária. 4 ed. São Paulo: Editora Manole. 1992.
3- WOLDEHIWET, Z.; RISTIC, M. Rickettsial and chlamydial diseases of domestic animals, 1ed, Oxford, Pergamon Press, 1993, 427 p.
4- WOODY, B.J.; HOSKINS, J.D. Erlichial Diseases of dogs. Veterinary Clinics of North America: Small Animal Practice. n.1. v. 21, p.129-140. Janeiro1991.
5- QUIBASA: Dados do Departamento de Pesquisa e Desenvolvimento.

QUALITY GUARANTEE

Before being released for consumption, all Quibasa reagents are tested by the Quality Control Department. The quality of the reagents is assured until the validity date mentioned on the packaged, as long as they are stored and transported under the appropriate conditions.
 **QUIBASA QUÍMICA BÁSICA Ltda**
Rua Teles de Menezes, 92 - Santa Branca
CEP 31565-130 - Belo Horizonte - MG - Brasil
Tel.: (31) 3439.5454
E-mail: bioclin@bioclin.com.br
CNPJ: 19.400.787/0001-07 - Industria Brasileira

CONSUMER SERVICE

Customer Support Service
Tel.: 0800 0315454
E-mail: sac@bioclin.com.br

Product licensed in the Ministry of Agriculture, Livestock and Supply since 31/01/2020 with the number 10.304/2020.

Technical manager: Dra. Camila Eckstein (CRMV/MG 20611)

Review: August/2021

UNIVERSAL SYMBOLOGY

CATALOG NUMBER

MANUFACTURED BY

BATCH CODE

CONTROL

DATE OF MANUFACTURE

POSITIVE CONTROL

USED BY
(last day of month)

NEGATIVE CONTROL

TEMPERATURE LIMITATION
(store at)

BIOLOGICAL RISK

CONTAINS SUFFICIENT
FOR <N> TESTS

INFLAMMABLE

CONSULT INSTRUCTIONS
FOR USE

CORROSIVE

IN VITRO DIAGNOSTIC DEVICE

POISON

EUROPEAN AUTHORIZED
REPRESENTATIVE

CE MARK

KEEP AWAY
FROM SUNLIGHT

DO NOT USE IF
PACKAGE IS
DAMAGED