Bioclin

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USAGE INSTRUCTIONS


FUNCTION

Method for the quantitative determination of Sodium in human samples of serum.. Enzymatic test, for *in vitro* Diagnostics use only.

PRINCIPLE OF ACTION

Methodology: Enzymatic

The Sodium is enzymatic determinate through the activity of the β -Galactosidase, depending of the Sodium with the ONPG as substrate. The absorbance at 405 nm of the product o-Nitrophenyl is proportional to the Sodium concentration.

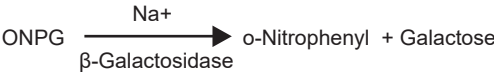ONPG = o-Nitrophenyl - β -D-Galactopyranose

REAGENTS

Number 1 - Buffer - Store between 2 and 8°C. Contains: β -D-Galactosidase < 8 U/mL, buffer, chelator and preservative.

Number 2 - Enzymatic Reagent - Store between 2 and 8°C. Contains: O-Nitrophenyl- β -D-Galactopyranose > 0,5 mM, buffer and preservative.

Number 3: Standard Level 1 - Store between 2 and 8°C. Contains: Sodium Chloride (concentration varies to each lot), buffers, surfactant, stabilizers and preservatives.

Number 4: Standard Level 2 - Store between 2 and 8°C. Contains: Sodium Chloride (concentration varies to each lot), buffers, stabilizers and preservatives.

PRESENTATION

Presenta-tion	Reagent Nº 1	Reagent Nº 2	Reagent Nº 3	Reagent Nº 4
1	40 mL	20 mL	3 mL	3 mL
2	2 x 40 mL	2 x 20 mL	3 mL	3 mL
3	3 x 40 mL	3 x 20 mL	3 mL	3 mL
4	4 x 40 mL	4 x 20 mL	3 mL	3 mL
5	30 mL	10 mL	3 mL	3 mL
6	2 x 30 mL	2 x 10 mL	3 mL	3 mL
7	3 x 30 mL	3 x 10 mL	3 mL	3 mL
8	4 x 30 mL	4 x 10 mL	3 mL	3 mL
9	60 mL	20 mL	3 mL	3 mL
10	2 x 60 mL	2 x 20 mL	3 mL	3 mL
11	3 x 60 mL	3 x 20 mL	3 mL	3 mL
12	4 x 60 mL	4 x 20 mL	3 mL	3 mL

EQUIPMENTS AND OPERATIONAL INPUTS

Biochemistry automated equipment, Biocontrol N and Biocontrol P Bioclin. These itens are found at markets specialized on Clinical Analysis Laboratories.

TRANSPORTATION AND STORAGE CONDITIONS

The storage temperature should be between 2 to 8°C. The transport at temperatures up to 30°C should not exceed 5 days. Protect from light and avoid moisture. **Do not freeze.**

SPECIAL CARE

- 1- For professional in vitro diagnostic use only.**
- Strictly follow the methodology proposed to obtain exact results.
- Water used in material cleaning must to be recent and free of contaminants.
- Do not ingest. Avoid contact with eyes and skin.
- The samples must be treated as potentially infectious.
- We recommend applying the local, state and federal rules for environmental protection, so that disposal of reagents and biological material can be made in accordance with current legislation.
- To obtain information related to biosafety or in case of accidents with the product, consult the MSDS (Material Safety Data Sheet) available on the website www.bioclin.com.br or upon request by the SAC (Customer Advisory Service) of Quibasa.
- Do not use the product in case of damaged packaging.
- It is essential that the instruments and equipments used are properly calibrated and subjected to periodic maintenance.

SAMPLES

Serum. The analyte is stable for 4 days between 2 and 8°C and 14 days at -20°C.

PROCESS DESCRIPTION

The calibration stability of the Sodium Ions Automated kit installed on refrigerated equipment is at least 3 days. This stability may vary depending on the conditions of the test, equipment and environment. Therefore, it is suggested to follow the product performance using control serum.

REAGENT PREPARATION

The reagents are ready for use.

TECHNIQUE

Bioclin recommends, as control serum, Biocontrol N and P Bioclin kits.

The kit is indicated for biochemistry automated analyzers only. **Check the application sheet for this equipment on the website www.bioclin.com.br or through SAC 0800 031 5454.**

RESULTS

The results are expressed in mmol/L.

PROCEDURE LIMITATIONS

Test periodically and calibrate to detect possible changes in the response of the equipment.

INTERFERENCE

The test does not suffer interference from the substances with the top concentrations of:

NH₄Cl 1.5 mmol/L
KPi 2.0 mmol/L
CaCl₂ 7.5 mmol/L
KCl 10 mmol/L
CuCl₂ 0.5 mmol/L
ZnCl₂ 0.5 mmol/L
FeCl₃ 0.5 mmol/L
Glucose 5 mmol/L
Ascorbate 10 mmol/L
Bilirubin 40mg/dL
Bilirubin conjugated 40 mg/dL
Hemoglobin 500 mg/dL
Triglycerides 1000 mg/dL

The anticoagulant Citrate is made from Sodium Citrate, then its use may cause falsely high results for the dosage of Sodium ion.

INTERNAL QUALITY CONTROL

The Clinical Laboratory must have an internal quality control, where all procedures, rules, limits and tolerance to variations be clearly established. It is important to mention that all measurement systems present a analytical variety, and it must be monitor by the laboratory. Therefore, it is recommendable the use of controls, allowing the precision and accuracy of the dosages.

TRACEABILITY

The kit's standards are traceable to the reference material NIST (National Institute of Standards and Technology) SRM 919B.

REFERENCE VALUES

The reference values for this method, were obtained through the determination of sodium in healthy populations of male and female.

Serum	136 - 146 mmol/L
	313 - 336 mg/dL

These values should be used as guidance, and each laboratory should establish its range of reference values, according to the population served. The results provided by this kit must be interpreted by the medical professional responsible, not being the only criterion for the determination of diagnosis and/or treatment of the patient.

PRODUCT PERFORMANCE
ACCURACY

Comparison of Methods

The Sodium Ions Automated was compared with another method commercially available to measure sodium. 42 analysis were conducted and the results were evaluated. The linear equation obtained was $Y = 0.9656X + 2.7094$, with a correlation coefficient 0.9805. With these results we can conclude that the kit has good methodological specificity.

PRECISION

Repeatability

The repeatability was calculated from 40 successive determinations, using 3 samples with different concentrations, obtaining the following results:

	Sample 1	Sample 2	Sample 3
Average Concentration (mmol/L)	125.23	141.33	160.33
Standard Deviation (mmol/L)	1.14	2.39	1.53
Coefficient of Variation (%)	0.91	1.69	0.95

Reproducibility

The reproducibility was calculated from 40 successive determinations for 3 consecutive days, using 3 samples with different concentrations, obtaining the following results:

	Sample 1	Sample 2	Sample 3
Average Concentration (mmol/L)	125.95	141.26	161.18
Standard Deviation (mmol/L)	0.63	0.14	0.76
Coefficient of Variation (%)	0.50	0.10	0.47

SENSITIVITY

The sensitivity was calculated from 40 determinations of a sample free of Sodium. The average found was -0.17 mmol/L with standard deviation of 0.99 mmol/L. The sensitivity, which indicates the method detection limit, corresponds the average plus 3 times the standard deviation and is equal to 2.79 mmol/L.

LINEARITY

The Sodium Ions Automated kit is linear between 80 and 180 mmol/L.

DIAGNOSTIC SIGNIFICANCE

In a healthy subject, the level of the Sodium extra cellular is regulated between 136-146 mmol/L (313-336 mg/dL). Small deviations of the regular level may cause great consequences in health. Sodium dosages are used in the diagnostic and treatment of aldosteronism, diabetes insipidus, adrenal hypertension, Addison disease, dehydration, inadequate secretion of the antidiuretic hormone, or other disease related to the electrolytic disequilibrium, metabolic and cardiovascular disorders. Therefore, the monitoring of the Sodium concentration in serum is important in both routine and emergency checking.

BIBLIOGRAPHIC REFERENCES

1. Berry, M. N. et al., (1988) Clin. Chem. 34,2295.
2. Tietz, N. W. (1983) Clinical guide to Laboratory Tests, p. 384 W.B. Saunders Co., Philadelphia.
3. QUIBASA: Dados do Departamento de Pesquisa e Desenvolvimento.

QUALITY ASSURANCE

Before being released for consumption, all **Bioclin** reagents are tested by the Department of Quality Control. The quality of reagents is assured until expiration date stated on the presentation packaging, when stored and transported under appropriate conditions.

QUIBASA QUÍMICA BÁSICA Ltda

Rua Teles de Menezes, 92 - Santa Branca
CEP 31565-130 - Belo Horizonte - MG - Brasil
Phone.: +55 31 3439.5454 | E-mail: bioclin@bioclin.com.br
CNPJ: 19.400.787/0001-07 - Made in Brazil

CUSTOMER SERVICE

Customer Advisory Service
Phone.: 0800 0315454 | E-mail: sac@bioclin.com.br

ANVISA registration for Sodium Ions Automated kit: 10269360211

Review: June/2022

UNIVERSAL SYMBOLOGY

CATALOG NUMBER

LOT NUMBER

MANUFACTURING DATE

VALIDITY DATE
(last day of the month)

TEMPERATURE LIMIT
(store)

CONTENT IS SUFFICIENT
FOR <N> TEST

SEE INSTRUCTIONS
FOR USE

IN VITRO DIAGNOSTIC
PRODUCT

KEEP AWAY
FROM SUNLIGHT

DO NOT REUSE

CAUTION

MADE BY

CONTROL

POSITIVE CONTROL

NEGATIVE CONTROL

BIOLOGICAL RISK

FLAMMABLE

CORROSIVE

TOXIC

DO NOT USE IF
PACKAGE IS
DAMAGED

PRODUCT
STERILIZED

DANGER